

AGS-Medical

Ph: 505-550-6501 & 505-565-5102

Fx: 505-814-5778, Em: sales@agsmedical.com

Web: <http://www.agsmedical.com>

WARNING

1. Please add water before using.
2. Away from flames and smoke.
3. When cleaning or repairing, you should cut off power supply.
4. Do not wet or spray liquids on the host because it would damage the equipment and cause electrical shock.



CE 0123



Eos-Ice Plus

The Golden Standard for Laser Hair Removal Platform--755nm, 808nm & 1064nm into a Single Innovative Handpiece



Eos-Ice Plus

More and more people desire to have permanent hair reduction. They want a treatment that is effective, fast and comfortable.

To meet these needs, a leader in laser industry, developed the Eos-Ice Plus.

Eos Ice Plus Platform offers the synergistic benefits of the 3 most effective wavelengths for hair removal, each targeting different structures within the hair follicle. The 3 main anatomical targets include the Bulge, Bulb and Papilla.

In order to provide you with optimal treatment possibilities, the vertical Eos-ice Plus comes with:

The "Gold standard" 808nm diode laser
Unique built-in USA semi-conductor cooling
Up to 20Hz repetition
USA military bar
Italy water pump
Laser power source TUV approved
Intuitive software with preset parameters and professional mode

The unique features of the Eos-Ice Plus, enable you to provide your patients with:

- ☑ Short treatment time
- ☑ Maximal treatment comfort
- ☑ Optimal results
- ☑ Skin type (I-VI) & hair color and texture adjustable

So You Can Grow with

Benefit from a rapid return on investment
Increase your patient satisfaction
Build up your clinic prestige
No topical anesthetic required
Reduce treatment times up to 45%
Differentiate your practice
Grow your profits

Eos - Ice Featuring Groundbreaking Dual Handle

755nm, 808nm & 1064nm into a Single
Innovative Handpiece.



Eos-Ice Plus



Technology behind Eos-Ice Plus

Celebrating 15 years of award-winning success, we are proud to present Eos Ice Plus Platform, a limited edition platform for 2017 featuring trio clustered diode technology. The Platform edition combines four laser wavelengths into a single innovative handpiece, simultaneously targeting different tissue depths within the hair follicle.

By combining the absorption and penetration levels of four different wavelengths, along with the treatment coverage, comfort and low maintenance, all kinds of patients can be treated without limitation of phototype, hair type or time of year with the maximum effectiveness and security.

Why choose Eos-Ice Plus ?

Premium USA Military Laser Generator

Alex 755nm: recommended for treating pale skins (phototypes I-III) with fine hair. Indicated for removing the residual hair of the last sessions. With more superficial penetration, the 755nm wavelength targets the Bulge of the hair follicle and is especially effective for superficially embedded hair in areas such as the eyebrows and upper lip.

Diode 808nm: recommended for treating all skin phototypes, particularly patients with a great density of hair and first sessions. Its deep penetration capabilities target the Bulge and Bulb of the hair follicle while moderate tissue depth penetration makes it ideal for treating the arms, legs, cheeks and beard.



Yag 1064nm: indicated for treating dark phototypes (III-IV tanned, V and VI). The 1064nm offers the deepest penetration of the hair follicle, allowing it to target the Bulb and Papilla, as well as treat deeply embedded hair in areas such as the scalp, arm pits and pubic areas.

Combined wavelengths: perfect balance to treat all phototypes (I-VI), we have the absorption efficiency of the 808nm wavelength, the power energy of the 755nm and the safety of the 1064nm wavelength. Excellent to arrive to all hair depths and dark skins.

Cooler Handle: Before and after treatment, cooler partner handle could be applied. Temperature could be adjusted on a handle. This handle can bring you new hair removal experience, easing the pain, enhancing the effect.

- ✓ Exclusive Patent
- ✓ Skin Preparation
- ✓ Skin Smoothing
- ✓ Laser Partner



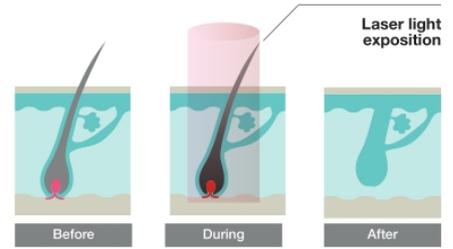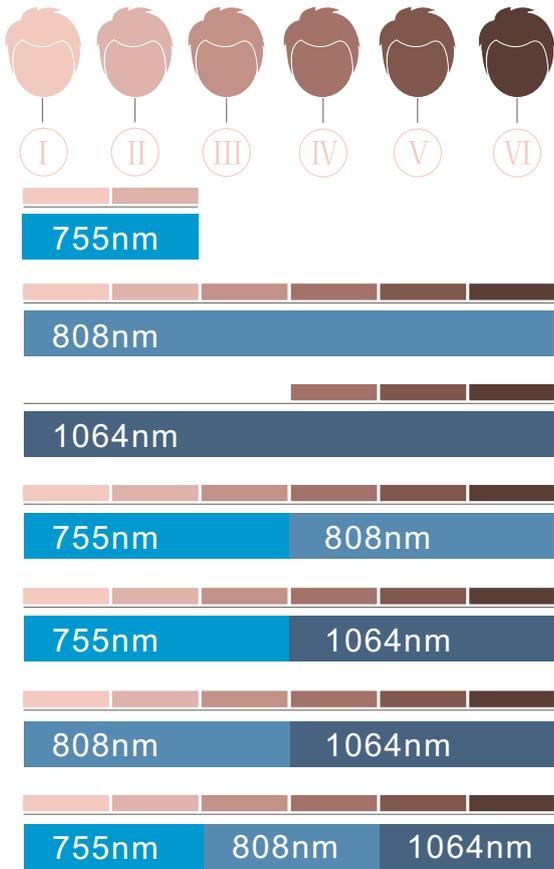

✔ Virtually Painless Treatment with Max 20HZ Speed

808nm diode laser adopts uniform energy that constantly penetrates deep into dermis and coagulate hair follicle with maximum 10 pulses in one second. Fast mode is specially designed to coagulate hair follicles to a proper damage point gradually with low fluence to ensure the painless and fast treatment.

✔ Friendly UI, Easy Menu Navigation

Presetting treatment proposals lead to easy navigation. Through many years' clinical study, diversified treatment proposals were concluded according to different skin types and skin areas. Besides, therapist could also freely set treatment parameters according to individual requirements.

Clinical Results

Clinically proven, over 95% of the cases have reported hair reduction of over 80% after 3 to 5 sessions of treatments. Hair become finer, lighter, fewer and grow slower. After 1 or 2 sessions, all cases reported a hair reduction of over 30% to 50%. The average energy density applied in most clinical cases is only 10 j/cm² to 20 j/cm² and sufficiently low energy avoids the potential risks.



SPECIFICATIONS

Eos-Ice Plus Platform	
Laser Type	Diode Laser
Laser Quantity	American military bars
Laser Wavelength	755 /808/1064/755+808/755+1064/808+1064/808+755+1064nm
The Light Guide Crystal	Sapphire
Spot Size	14x14mm
Repetition Frequency	1 ~ 20HZ
Pulse Width	10-400ms
Energy Density	1 ~ 120j/cm ²
Laser Cooling Temperature	-5°C - 5 °C
Power Supply	AC230V,50Hz / AC110V, 60Hz
Cooling System	Air + Water + Semiconductor + Sapphire
Laser Cooling Water Temperature	20 - 30°C
Power Consumption	2000VA
Rated Power	3000W
Fuse	Model: F250V L15A
Package Size	Machine: 109*62*48cm / Handpiece: 33*46*12cm
Net Weight	Machine: 49kg / Handpiece: 1.9kg

WARNING

1. Please add water before using.
2. Away from flames and smoke.
3. When cleaning or repairing, you should cut off power supply.
4. Do not wet or spray liquids on the host because it would damage the equipment and cause electrical shock.



CE 0123

